

Report_andrea_braze Feb 28 2024 12:38:19 GMT+1 OPENSIGNALS VERSION: Public Build 2023-10-23

Author: PLUX


Devices

	name OC:43:14:1C:2A:25			bitalino_rev	1000 Hz
channel label sensor	resolution	channel	label	sensor	resolution
Al EMGBITREV	10	П	INPUT	BTN	1
A2 A2 RAW	10	I 2	INPUT	Generic Input	1
A3 RAW	10	01	ОИТРИТ	LED	1
A4 A4 RAW	10	02	ОИТРИТ	BUZ	1
A5 A5 RAW	10	DAC	DAC	PWM	8
A6 A6 RAW	10				





Electromyographic analysis

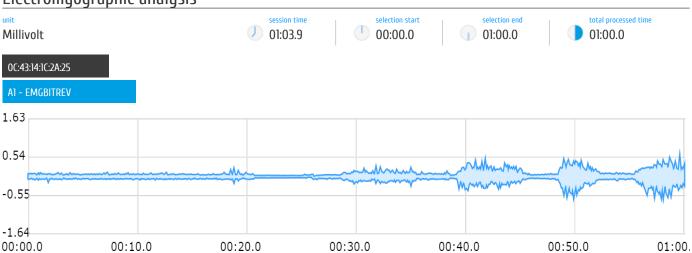

Processed channels

0C:43:14:1C:2A:25

A1

FILTER

applied bandpass - [10,100] Hz cutoff - Order: 6

DETECTION OPTIONS

applied to TKE: smooth window of 0.5 seconds

Muscle Activation results

Activations

find all muscle activations

0C:43:14:1C:	2A:25												
channel	digital	activation	start (s)	end (s)	duration (s)	latency (s)	max (mV)	min (mV)	ampl+ (mV)	ampl- (mV)	ampl+- (mV)	rms (mV)	area (mV^2)
1	NONE	1	29.493	30.109	0.616	29.493	0.13 mV	-0.19 mV	25.00	35.00	0.32 mV	0.04 mV	20 mV ²
	NONE	2	30.820	30.995	0.175	30.820	0.13 mV	-0.12 mV	24.00	24.00	0.25 mV	0.04 mV	5.28 mV ²
	NONE	3	31.674	31.814	0.140	31.674	0.2 mV	-0.16 mV	38.00	30.00	0.36 mV	0.05 mV	3.94 mV ²
	NONE	4	32.084	33.050	0.966	32.084	0.2 mV	-0.18 mV	38.00	34.00	0.38 mV	0.05 mV	32.58 mV ²
	NONE	5	34.841	35.168	0.327	34.841	0.19 mV	-0.22 mV	36.00	41.00	0.41 mV	0.04 mV	9.41 mV ²
	NONE	6	39.227	42.609	3.382	39.227	0.38 mV	-0.37 mV	71.00	70.00	0.75 mV	0.07 mV	172.4 mV ²
	NONE	7	42.644	45.168	2.524	42.644	0.31 mV	-0.31 mV	58.00	58.00	0.62 mV	0.06 mV	101.96 mV ²
	NONE	8	48.493	51.955	3.462	48.493	0.42 mV	-0.55 mV	79.00	103.00	0.96 mV	0.09 mV	223.18 mV ²
	NONE	9	55.620	59.957	4.337	55.620	0.47 mV	-0.55 mV	88.00	103.00	1.02 mV	0.1 mV	293.23 mV ²
Diff	-	-	-	-	-	-	-	-	-	-	-	-	-
Ratio	-	-	-	-	-	-	-	-	-	-	-	-	-

Glossary



^{1.} MVC: Maximum Voluntary Contraction.

^{2.} TKE: Teager-Kaiser Energy operator.



Signal statistics analysis total processed time 0.00:00 01:03.9 01:00.0 01:00.0 0C:43:14:1C:2A:25 A2 - EMGBITREV 1.63 0.54 -0.55 -1.64 00:00.0 00:10.0 00:20.0 00:30.0 00:40.0 00:50.0 01:00. A3 - RAW 1023.00 682.00 341.00 0.00 0.00:00 00:10.0 00:20.0 00:30.0 00:40.0 00:50.0 01:00. A4 - RAW 1023.00 682.00 341.00 0.00 00:00.0 00:10.0 00:20.0 00:30.0 00:40.0 00:50.0 01:00. A5 - RAW 1023.00 682.00 341.00

00:30.0

00:40.0



00:50.0


01:00.

00:10.0


00:20.0

0.00



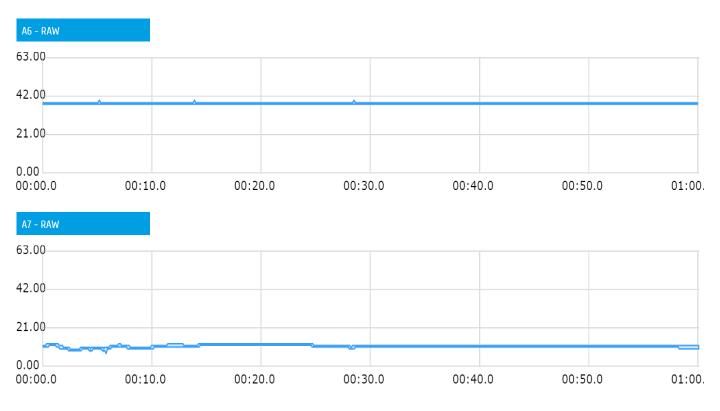

Processed channels

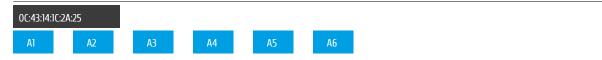

Signal Statistics

0C:43:14:1C:2A:25										
CHANNEL	MIN	MAX	MEAN	MEDIAN	STD	PEAK2PEAK	RMS	RSS	UNIT	
1	-0.55	0.47	-0.01	-0.02	0.04	1.02	0.05	11.11	mV	
2	504.00	518.00	510.69	511.00	1.76	14.00	510.69	125092.95		
3	0.00	0.00	0.00	0.00	0.00	0.00	0.00	0.00		
4	0.00	1019.00	520.13	511.00	76.84	1019.00	525.78	128788.47		
5	38.00	39.00	38.00	38.00	0.01	1.00	38.00	9308.10		
6	8.00	12.00	11.06	11.00	0.63	4.00	11.07	2712.58		

